

Operating Manual Ionizer

KERN YBI-01

Version 1.1
03/2013
GB



YBI-01-BA-e-1311



KERN YBI-01

Version 1.1 03/2013

Operating Manual Ionizer

Table of Contents

1	GENERAL	4
2	BASIC SAFETY PRECAUTIONS	4
3	TECHNICAL DATA	7
4	APPLIANCE OVERVIEW	8
5	UNPACKING, SETUP AND COMMISSIONING	9
	Unpacking	9
	Placing.....	9
	Mains connection	11
6	COMMISSIONING	11
	Turning on the ionizer.....	11
	Turning on the fan.....	11
	Applications:.....	12
7	SERVICE, MAINTENANCE, DISPOSAL	13
	Cleaning.....	13
	Maintenance Service.....	13
	Disposal.....	13
8	DECLARATION OF CONFORMITY	14

1 GENERAL

The ionizer has conductive peaks supplied by high voltage which, based on corona discharge, generate positively and negatively charged ions in the immediate vicinity. These are attracted by the electrostatic charge of goods to be weighed and, thus, neutralise the interfering electrostatic charge. This also does away with the forces falsifying the weighing (such as falsified weighing result, weighing value drifted).

2 Basic Safety Precautions

WARNING



The use of ionizer is only intended in combination with electronic weighing balances. Do not use for any other purposes.



Never operate the ionizer in explosive environment. The serial version is not explosion protected.



Protect the ionizer against high air humidity / temperature, steams and dust;

Take care to select a location free of water and oil

Do not expose the ionizer to strong humidity for extended periods. Non-permitted condensation (condensation of air humidity on the appliance) may occur if a cold appliance is taken to a considerably warmer environment. In this case, acclimatize the disconnected ionizer for ca. 2 hours at room temperature.



Whilst the ionizer is switched on, do not touch the ion source; see sticker on the left side.



In case of smoke development, smell of fire, strong heating-up of the ionizer or when the red LED starts glowing, turn off the ionizer immediately and disconnect it from the mains.



If water or other foreign matter enters the ionizer, turn off the master switch immediately and disconnect it from the mains.



On account of the high voltage technology, handle the ion source and exits with care.



Do not take apart or modify the ionizer.



Prevent damage caused by dropping, vibration or shock; see sticker on the left side.



Always use the genuine power pack. The stated voltage value must be the same as the local voltage.



Risk of injury! The peaks of the ion source are sharp and cutting.



The ionizer generates poisonous ozone; ensure sufficient ventilation.



For maintenance and repair work disconnect the ionizer from the mains.



Disconnect the ionizer from the mains during periods of idleness.

CAUTION



Maintain and clean the ionizer at regular intervals.

- Clean the ion source after 100 hours.
- Replace the ion source after 10 000 hours.



Starting up a damaged ionizer may result in a short circuit, fire or electric shock.



Starting up out of doors and inside vehicles is prohibited and will result in total loss of warranty.



The occurrence of electromagnetic fields may result in major display deviations (incorrect weighing results). Discharge sample at a sufficient distance from the weighing balance.



During normal operation the green LED will be glowing, in case of a breakdown the red LED.

If the red LED is glowing, turn off the ionizer at the master switch and turn it on again. If the red LED continues to glow, inform the manufacturer.



The ionizer may heat up slightly during operation.



Prior to ionising powdery goods to be weighed, turn off the fan.




Ionisation results are better when the fan is turned on and the period required for discharging the sample will be shortened.

3 Technical Data

Distance "sample ion source"	c. 5 – 50 cm
Discharge time	8 s/5 cm, 12 s/10 cm, 100 s/50 cm (with fan turned on)
Ozone concentration	0.04 ppm (2cm from ion source)
Ion source	Stainlees steel
Maximum air volume	0.06 cm ³ /min
Weight	Ionizer 110 g, frame 430 g
Ambient conditions	0- 40 °C, 35 – 80 % air humidity (non-condensing)
AC adapter (primary)	100 – 240V, 400 mA 50/60Hz
Rated electric power supply	DC 12V, 1. A
Pollution Degree	2
Overvoltage Category	Category II
Altitude	Up to 2000 m
Installation Site	device may only be used indoors

4 Appliance overview

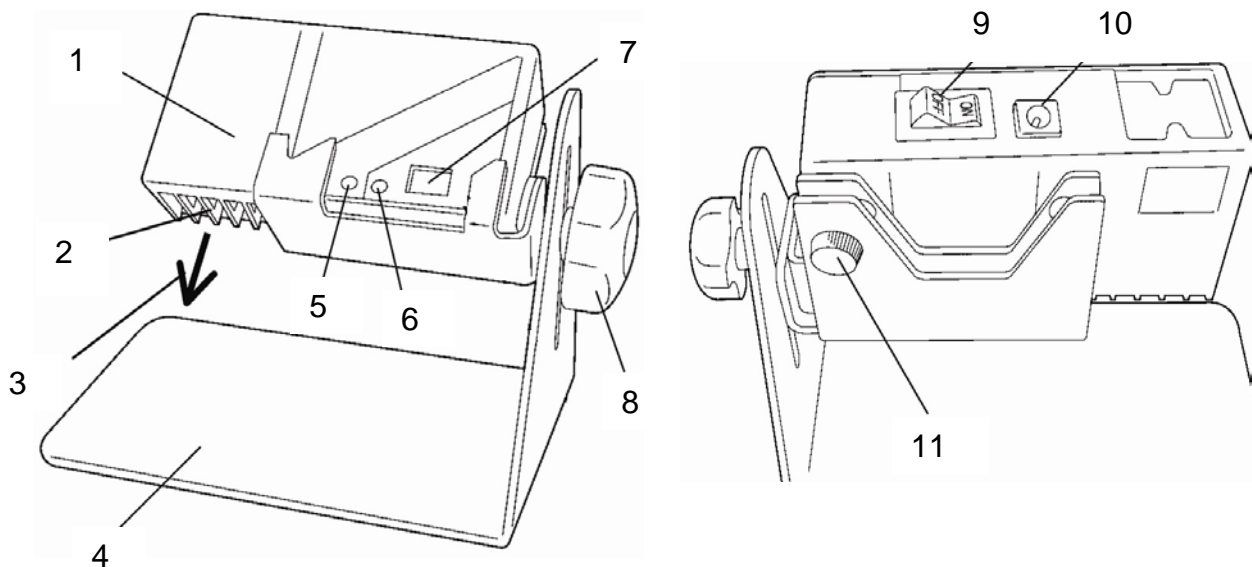


1. Ionizer
2. Ion source
3. Discharge direction
4. Frame
5. LED Ionizer operation (green = normal operation, red = breakdown)
6. LED fan operation (blue = on)
7. Fan on- and off-switch
8. Knob for setting height / angle of inclination
9. Main switch
10. Mains adapter connection
11. Knob for permanent fixing of ionizer

5 Unpacking, Setup and Commissioning

Unpacking

Remove device carefully from packaging, remove plastic wrapping and place device at the planned work place.

Scope of delivery

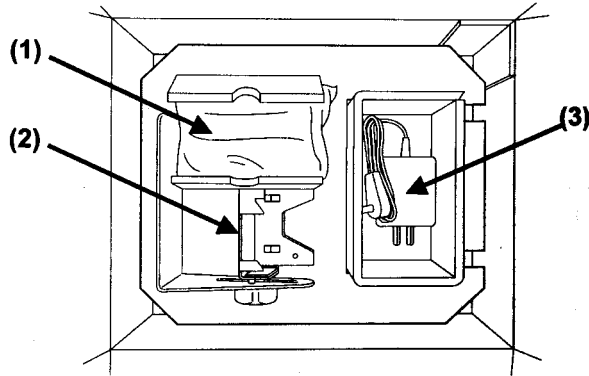


- (1) Ionizer
- (2) Frame
- (3) Mains adapter
- (4) Operating instructions

Placing


The ionizer may be used as a desktop or hand-held device.

Using it on a desktop

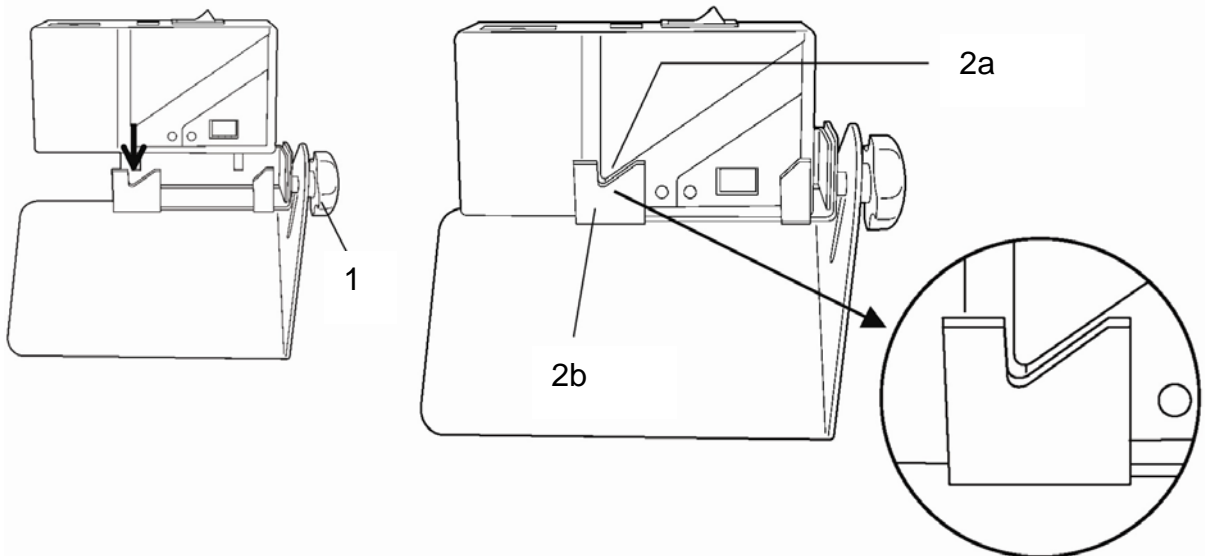

Position knob (1) and tighten it by turning.
Insert the ionizer into the V-shaped recess (2a → 2b),
ensuring a good fit


If the device is to be used permanently on a desktop, fix it with a set screw (3).

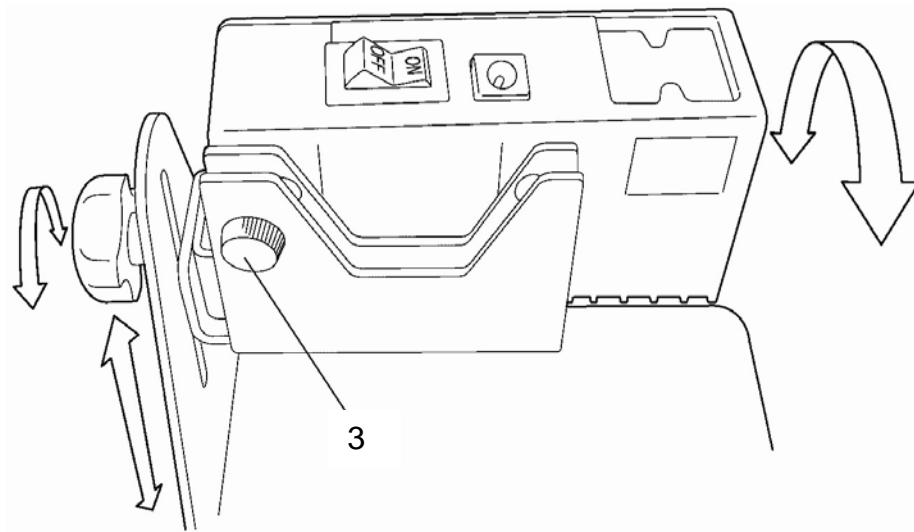

Use the knob on the frame to set the height and angle of inclination.

Using it as a hand-held device

To use it as a hand-held device, remove the ionizer from the frame as illustrated.

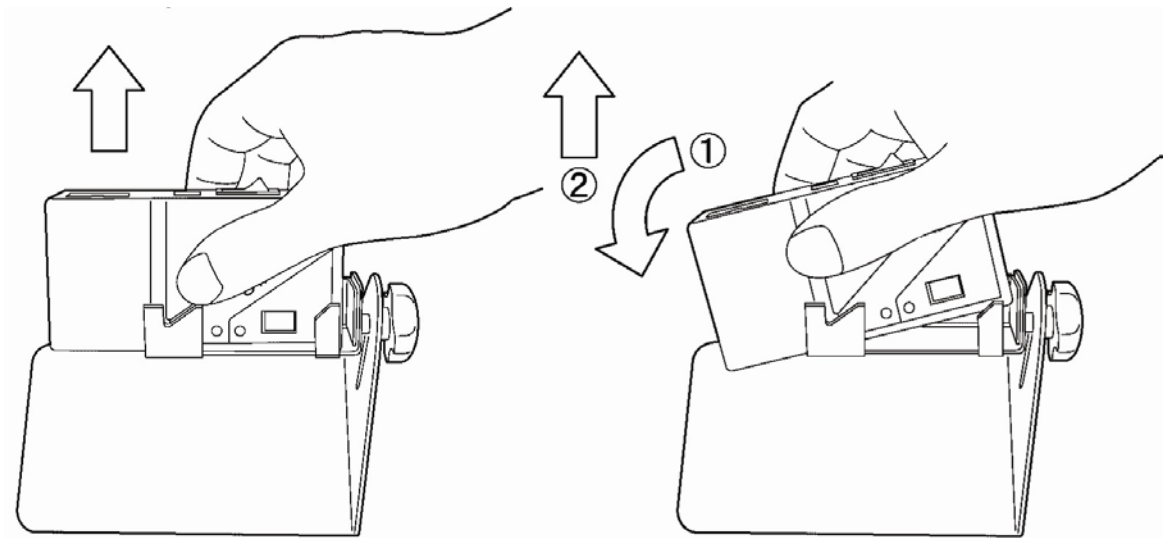

Mains connection

Power is supplied via the external mains adapter. The stated voltage value must be the same as the local voltage.

Only use original KERN mains adapters. Using other makes requires consent by KERN.

Turn off the ionizer (master switch "OFF") before connecting it to the mains adapter.



6 Commissioning

Turning on the ionizer

- ⇒ Set master switch to "ON".
The green LED starts glowing.

During normal operation the green LED will be glowing, in case of a breakdown the red LED.

If the red LED is glowing, turn off the ionizer at the master switch and turn it on again. If the red LED continues to glow, disconnect the device from the mains and contact the manufacturer.

After use turn off with master switch "OFF".



Turning on the fan

Ionisation results are better when the fan is turned on and the period required for discharging the sample will be shortened.

Information:

Turn off fan for powdery, easily whirled up samples.

- ⇒ Set switch of fan to "ON".
The blue LED starts glowing.

Ensure sufficient ventilation even when the fan is turned on.



Applications:

The ionizer is designed to be used in combination with electronic weighing balances!



Use as desktop device.

⇒ Discharge of solids or weighing containers.

Ionisation results are better when the fan is turned on and the period required for discharging the sample will be shortened.



Use as desktop device.

⇒ Discharge of powdery samples. Discharge prevents whirling up, problems with toxic samples.


Turn off fan for easily whirled up samples.



Use as hand-held device inside weighing chamber.

⇒ Discharge of goods to be weighed, glass wind screen or weighing containers.

If your hand is inside the weighing chamber it will generate heat dissipation. Wait 2-3 minutes before reading the weighing value.

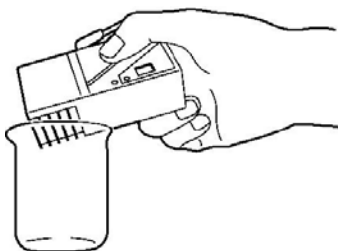

Use as hand-held device outside the weighing chamber.

⇒ Discharge of beakers etc.

Powdery samples are deposited at the inner wall of the glass when being poured into the beaker. This is prevented by discharging the beaker.

7 Service, maintenance, disposal

Cleaning

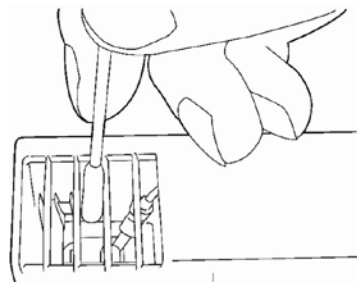


Fig. 1:

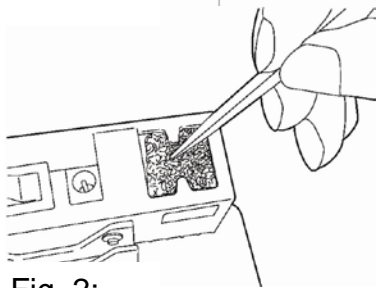


Fig. 2:

WARNING



Disconnect the device from the mains before cleaning / change of filter.



Do not take apart the ionizer.



Insert only completely dried filters.

Do not use aggressive cleaning agents (solvents etc.) instead clean with a soft cloth soaked in mild soap suds. Take care that the device is not penetrated by fluids and polish it with a dry soft cloth.

Remove loose sample residues / powder carefully with the help of a brush or hand-held vacuum cleaner.

For cleaning the ion source use a cotton stick moistened with alcohol; see fig. 1. Take care that the peaks are not bent.

For cleaning the filter use tweezers to remove it and clean it with water; see fig. 2. Wait until completely dry and reinsert.

Maintenance Service

The appliance may only be opened by trained service technicians who are authorized by KERN.

Before opening, disconnect from power supply.

Disposal

Disposal of packaging and appliance must be carried out by operator according to valid national or regional law of the location where the appliance is used.

8 Declaration of conformity



KERN & Sohn GmbH

D-72322 Balingen-Frommern
Postfach 4052
E-Mail: info@kern-sohn.de

Tel: 0049-[0]7433- 9933-0
Fax: 0049-[0]7433-9933-149
Internet: www.kern-sohn.de

Declaration of conformity

EG-Konformitätserklärung

EC- Déclaration de conformité

EC-Dichiarazione di conformità

EC- Declaração de conformidade

EC-Deklaracja zgodności

EC-Declaration of -Conformity

EC-Declaración de Conformidad

EC-Conformiteitverklaring

EC- Prohlášení o shode

ЕС-Заявление о соответствии

D	Konformitäts- erklärung	Wir erklären hiermit, dass das Produkt, auf das sich diese Erklärung bezieht, mit den nachstehenden Normen übereinstimmt.
GB	Declaration of conformity	We hereby declare that the product to which this declaration refers conforms with the following standards.
CZ	Prohlášení o shode	Tímto prohlašujeme, že výrobek, kterého se toto prohlášení týká, je v souladu s níže uvedenými normami.
E	Declaración de conformidad	Manifetamos en la presente que el producto al que se refiere esta declaración está de acuerdo con las normas siguientes
F	Déclaration de conformité	Nous déclarons avec cela responsabilité que le produit, auquel se rapporte la présente déclaration, est conforme aux normes citées ci-après.
I	Dichiarazione di conformità	Dichiariamo con ciò che il prodotto al quale la presente dichiarazione si riferisce è conforme alle norme di seguito citate.
NL	Conformiteit- verklaring	Wij verklaren hiermede dat het product, waarop deze verklaring betrekking heeft, met de hierna vermelde normen overeenstemt.
P	Declaração de conformidade	Declaramos por meio da presente que o produto no qual se refere esta declaração, corresponde às normas seguintes.
PL	Deklaracja zgodności	Niniejszym oświadczamy, że produkt, którego niniejsze oświadczenie dotyczy, jest zgodny z poniższymi normami.
RUS	Заявление о соответствии	Мы заявляем, что продукт, к которому относится данная декларация, соответствует перечисленным ниже нормам.

Ionizer: KERN YBI-01

EU Directive	Standards
2004/108/EC	EN 61326-1: 2006
2006/95/EC	EN 61010-1: 2010

Datum 27.03.2013

Date

Signatur

Signature

Ort der Ausstellung 72336 Balingen

Place of issue

Albert Sauter
KERN & Sohn GmbH
Geschäftsführer
Managing director

KERN & Sohn GmbH, Ziegelei 1, D-72336 Balingen, Tel. +49-[0]7433/9933-0
Fax +49-[0]7433/9933-149, E-Mail: info@kern-sohn.com, Internet: www.kern-sohn.com